

Service Information Bulletin

This Bulletin is for information only.

Release of new Sparklan Driver to support EAP

Task

Timing		Category
Next service	<input type="radio"/>	Apply at all sites
	<input checked="" type="radio"/>	Apply at affected sites as listed below
	<input type="radio"/>	Optional to improve functionality of product

Task Tracking

After completion of your task the following entry in your Service Report is required:

DD+DIS009.15E

*

* Insert the document number into the field "Comment" (SMS form).

Purpose of this document:

- This document informs the field that a new Sparklan driver is available.
- The new Sparklan driver version will support Extensible Authentication Protocol (EAP) authentication in WiFi networks.
- This document provides download information for the driver.

Affected serial number(s) / batch:

1 Introduction/purpose

Symptom EAP authentication for WiFi networks is not possible at all hospital sites where DX-D 100 / 100 Wireless is used.

Cause EAP authentication is not supported by the current Sparklan driver.

Solution Install the new Sparklan driver as described in section 3.



NOTE:

The new driver only supports IEEE 802.11g and IEEE 802.11n protocols (2.4GHz & 5GHz).

It does not support IEEE 802.11b.

2 Prerequisites



IMPORTANT:

Configuration of the EAP authentication protocol depends on the site infrastructure. The local IT manager needs to be consulted before installation of the firmware.



SOFTWARE:

Download the Mediatek wireless utility from the Agfa HealthCare Library. It includes the Sparklan driver RT2870 V 5.1.17.0:

Direct Radiography > DX-D 100 Wireless > DX-D 100 Wireless Generic > Software



TOOLS:

Standard service engineers' tool kit

Flash storage device for storing the Sparklan Driver RT2870 V 5.1.17.0

3 Instructions



REQUIRED TIME:

Approx. 30 minutes:

- Approx. 10 minutes installation of the driver
- Approx. 20 minutes verification

Install the driver for the SparkLAN WUBR-507N WiFi network card as described in the following steps:

- (1) Stop NX.
- (2) Insert the Flash drive medium at USB port of the built-in NX PC. The port is situated at the Peripheral Connections Panel. See red mark in figure 1.
- (3) Select the Mediatek wireless utility from the Flash drive.

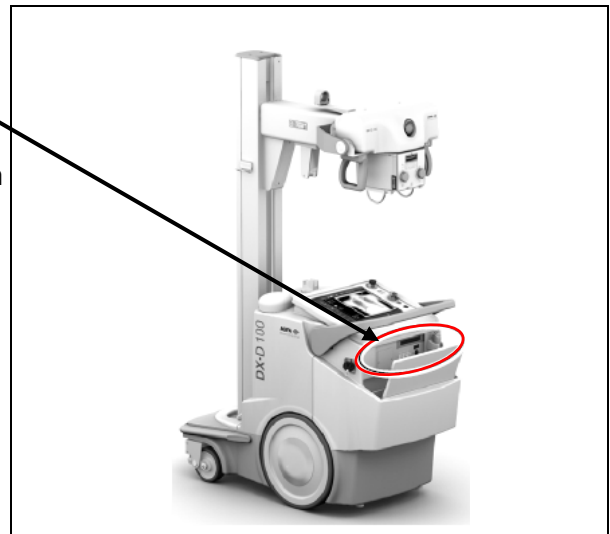


Figure 1

- (4) Double-click the exe inside:
IS_Setup_ICS_081414_1.5.38.101.exe
The screen shown in figure 1 appears:

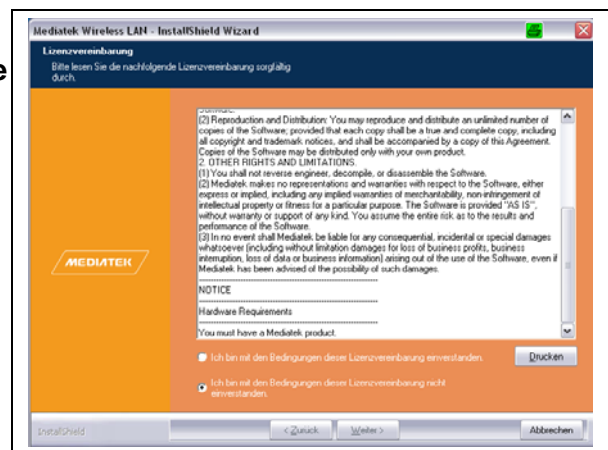


Figure 2

- (5) Accept the conditions of the license agreement. See blue frame in figure on the right.

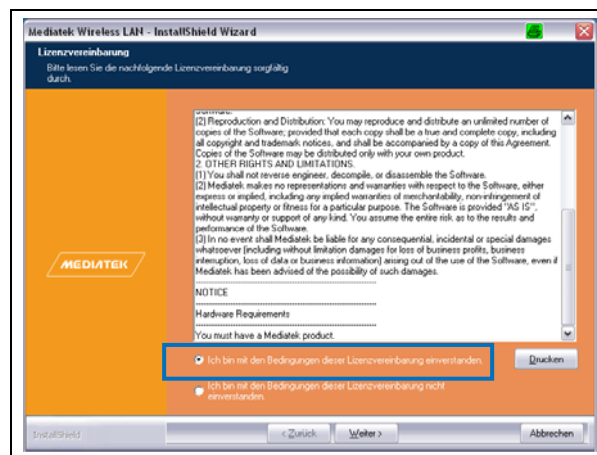


Figure 3

- (6) Click **Next**
(7) Select **install driver only**. See blue frame in figure on the right.

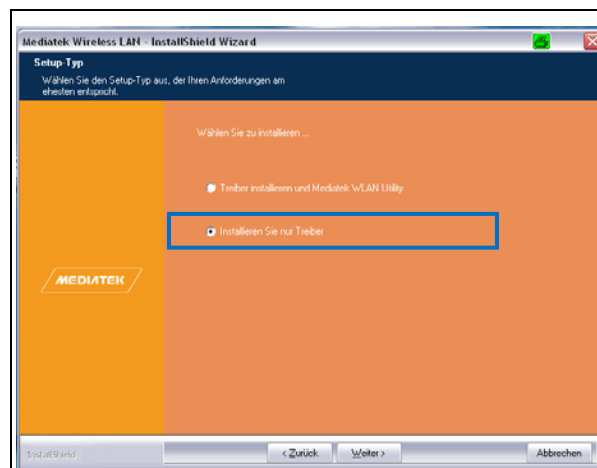


Figure 4

- (8) Click **Next**.
(9) Click **Install**.
(10) Install the driver.
(11) Reboot NX.

4 Verification

Check if the IEEE802.11g and the IEEE802.11n protocols are supported after installation of the new driver.

Check if the connection to the hospital network is working.

5 Keywords

WiFi, EAP, IEEE802.11n protocol, IEEE802.11g protocol, Sparklan, driver

6 Version history

Version	Change	Date
1.0	Initial version	02-2015