

DR 100e	Type 6012	No. 51
DR 100s	Type 6013	No. 50
DR 400	Type 5520	No. 169
DR 600	Type 5530	No. 215
DR 800	Type 6010	No. 60
DR Retrofit	Type 5400	No. 134
DX-D 100	Type 5410	No. 141
DX-D 100 Wireless	Type 5411	No. 190
DX-D 300	Type 8207	No. 175
DX-D 400	Type 5420	No. 129
DX-D 600	Type 5430	No. 196
Valory	Type 5540	No. 13

Document ID: 80106998

Service Information Bulletin

This Bulletin is for information only.

Service Quality Test Tool version 2.3 released

Task

Timing		Category	
<input checked="" type="radio"/>	Next service	<input type="radio"/>	Apply at all sites
		<input checked="" type="radio"/>	Apply at affected sites as listed below
		<input type="radio"/>	Optional to improve functionality of product

Task Tracking

After completion of your task the following entry in your Service Report is required:

*

* Insert the document number into the field "Comment" (SMS form).

Purpose of this document:

This document announces the release of Service Quality Test (SQT) Tool version 2.3.

Affected systems:

The Service Quality Test Tool is used for all DR systems listed on the front page of this document.

1 Introduction

Version 2.3 of the Service Quality Test Tool (SQT) is available for download.

It replaces the previous version 2.2.

SQT 2.3 is updated with XD 17 and XD+17 detectors used as fixed detectors in a Bucky. It is recommended to use only the SQT 2.3 (or later version) from now on.



NOTE:

For more information to the SQT Tool see also:

- Service Bulletin *Service ATPs replaced by unified Service Quality Test Tool*, Document ID [54799858](#).

2 Prerequisites



SOFTWARE:

- Service Quality Test Tool 2.3 ZIP File, Document ID [55022710](#), comprising:
 - SQT_2.3.0000.exe incl. test instructions* (Online help)
 - Using the SQT.PDF (short instructions)
- Excel ≥ 2007
- Adobe Acrobat Reader or Adobe Acrobat.
Other PDF readers in general also work, but they will not open the SQT instructions with the specific page for the test.

Download the SQT 2.3 from the Agfa Medimg Library or directly from the AgfaBox. For details refer to the Service Information Bulletin *AgfaBox replaces Fileshare for software downloads*, Document ID [71987422](#).

* SQT instructions can also be downloaded separately with Document ID [54513844](#).



IMPORTANT:

The installation and service of the product(s) described herein is to be performed by qualified personnel who are employed by Agfa NV or one of its affiliates or who are otherwise authorized by Agfa NV or one of its affiliates to provide such services.



IMPORTANT:

Strictly observe all safety directions within the *Generic Safety Directions*, Document ID [11849633](#).

3 Instructions



REQUIRED TIME:

- 5 minutes: Download of SQT 2.3
- 1 minute: SQT 2.3 installation
- 5 minutes: Installation of .NET Framework 4.8 if not yet available

- (1) Download and install the SQT 2.3.
SQT 2.3 automatically uninstalls previous SQT versions and preserves the available SQT site and system database.

The SQT is installed in directory C:\Agfa\HealthCare\SQT.



IMPORTANT:

- The list of “Incomplete Test Runs” is not preserved during SQT upgrade.
- When installing the same or previous SQT version, the “site and system database” (function “Select from History”) is deleted.
Reason: When the same or previous version is installed, it is assumed that it was intended to repair the SQT.


4 Verification



IMPORTANT:

After a service intervention, make sure that the device is operational as intended, before you hand it over to the customer.
Follow the instructions in this section precisely.

- (1) Start the SQT 2.3 using the shortcut in the start menu.

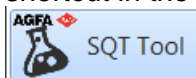

- (2) Confirm that the SQT 2.3 start page appears.

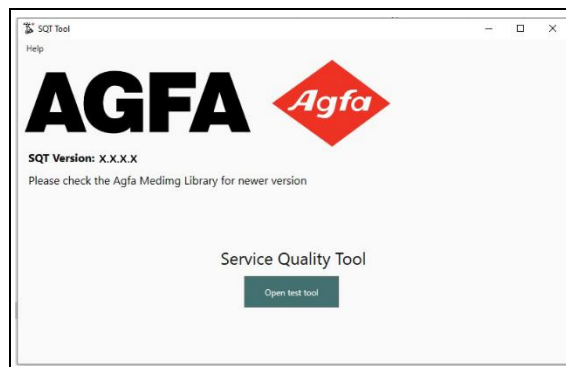


Figure 1

- (3) If error message **This application requiresNETFramework,Version=v4.8** appears, click **Yes** and follow the instructions to install it.
This error should only appear for new SQT installations.

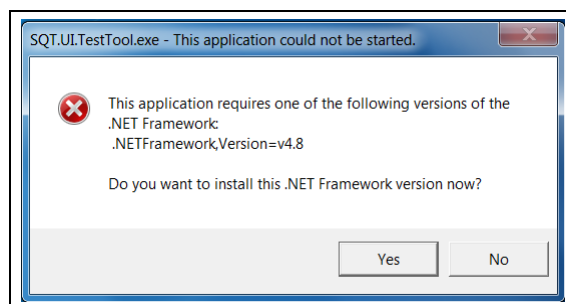


Figure 2

5 Keywords

acceptance, ATP

6 Version history

Version	Change	Date
1.0	Initial Version	08-2022