

## Service Information Bulletin

This Bulletin is for information only.

### Revised Information:

Information added in case of an exchange  
from Crystal-lead battery to Lead-acid batteries

## Re-introduction of the Lead-acid batteries as of January 2021 in production

### Task

Timing		Category
<input type="radio"/>	Next service	<input checked="" type="radio"/> Apply at all sites
		<input type="radio"/> Apply at affected sites as listed below
		<input type="radio"/> Optional to improve functionality of product

### Task Tracking

After completion of your task the following entry in your Service Report is required:

75158226

\*

\* Insert the document number into the field "Comment" (SMS form).

Purpose of this document:

This document contains all necessary information about:

- DX-D 100 and DX-D 100 Wireless will be produced from January 2021 on again with Lead-acid batteries as the battery manufacturer stopped production for the Crystal-lead batteries.
- Crystal-lead batteries can be replaced by Crystal-lead batteries only, so spare parts had been put on stock to cover the next ten years.

Affected sites:

Lead-acid batteries are re-introduced January 2021 in production:

Use, dissemination, distribution or reproduction of this bulletin by unauthorized personnel is not permitted and may be unlawful.  
The controlled version of this document is available from the Agfa Medimg Library. Any printed copy of this document is uncontrolled.

Copyright © 2021 Agfa NV | Document ID: 75158226

## 1 Introduction/purpose

The Battery manufacturer stopped the production of the Crystal-lead batteries.



**NOTE:**

This will be effective for DX-D 100 and DX-D 100 Wireless as of January 2021.



**IMPORTANT:**

Availability of spare parts of Crystal-lead battery powered systems is secured for 10 years after the date of manufacturing of the last system with crystal lead batteries.

From January 2021 the Lead-acid batteries are introduced as successor in production.

## 2 Prerequisites



**SPARE PARTS:**

For the related spare parts refer to:

DX-D 100 / DX-D 100 Wireless - Mobile X-Ray Unit - Spare Parts List,  
Document ID [30872155](#)

In case of an exchange from Crystal-lead battery to Lead-acid batteries order:  
SC+A521178-01 CONVERSION KIT TO LEAD-ACID BATTERIES



**NOTE:**

The CONVERSION KIT is on “manual distribution”, it requires an approval from HQ to get it delivered.

The CONVERSION KIT contains the following parts:

A. Supplement Motor Battery Tray	2 Units
B. Battery Tray W. 15Ah Battery	5 Units
C. Battery Monitor SAT-A3580-12	1 Unit
D. Tie-Wraps	7 Units
E. Warning Label	1 Unit
F. Battery Connectors Labels	1 Unit
G. Screw 4X8	4 Units
H. Screw 4x16	4 Units

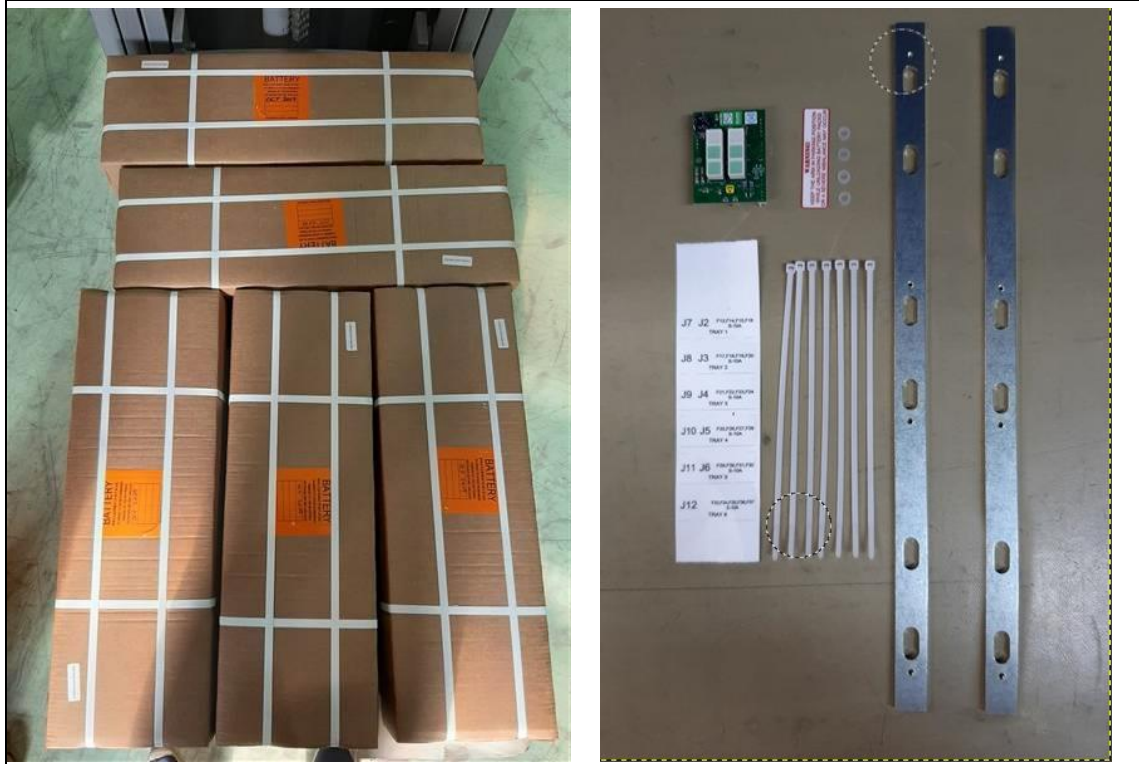


Figure 1: Scope of delivery



**IMPORTANT:**

The installation and service of the product(s) described herein is to be performed by qualified personnel who are employed by Agfa NV or one of its affiliates or who are otherwise authorized by Agfa NV or one of its affiliates to provide such services.



**IMPORTANT:**

Strictly observe all safety directions within the "Generic Safety Directions", Document ID [11849633](#).

### 3 Instructions



REQUIRED TIME:

8 hours

- (1) Follow the instructions as described in:



SIN 19-07-08 Rev.C  
UPGRADE ADAPTATI

### 4 Verification

- (1) Perform a DR workflow to verify functionality, charge level and performance of the new batteries.
- (2) Inform customer to charge the equipment for at least for 10h.

### 5 Keywords

Generator, battery, crystal lead batteries, lead-acid batteries

### 6 Version history

Version	Change	Date
1.0	Initial Version	12-2020
1.1	Revised Information: Information added in case of an exchange from Crystal-lead battery to Lead-acid batteries	04-2021