

Service Bulletin

New WLAN Adapter Available for Solving Performance and Connection Reliability Issues

Task			
Timing		Category	Scope
Next service as agreed with customer	<input type="radio"/>	Apply at all sites	PowerHelp complaint / HQ issue: HQ_1110270002
	<input checked="" type="radio"/>	Apply at affected sites as listed below	
	<input type="radio"/>	Optional to improve functionality of product	

Task Tracking

After completion of your task the following entry in your Service Report is required:

DD+DIS179.12E

*

* Insert the document number into the field "Comment" (SMS form).

Purpose of this document:

This document informs about the availability of a new WLAN Adapter for the DX-D 100 solution.

This WLAN Adapter will support the IEEE802.11n protocol to ensure operation in all hospital environments.

The upgrade to Windows 7 improves log-off and battery conservation functionalities

Affected serial number(s) / batch: This Service Bulletin is only applicable for XP-based DX-D 100 systems.

The signatures on the approval page indicate the solutions described in this Service Bulletin have been reviewed and are NOT reportable because no actions are taken to reduce a "Risk to Health" according to our risk assessment process.

NOTE:



This Bulletin replaces the earlier Service Bulletin No. 11 - DD+DIS086.12E, External USB Adapter available for Solving Performance and Connection Reliability Issues.

Use, dissemination, distribution or reproduction of this bulletin by unauthorized personnel is not permitted and may be unlawful. The controlled version of this document resides on the HealthCare Library. Any printed copy of this document is uncontrolled.

Copyright © 2009 Agfa HealthCare N.V. | Document Node ID:37440705 | service_bulletin_template_v06

1 Introduction/purpose

Symptoms Several (connectivity) issues with the current WLAN Adapter on XP-based DX-D100 installations.

Cause The IEEE802.11n protocol was not supported by the temporary solution (see Service Bulletin No. 11 - DD+DIS086.12E).

Solution Install the Upgrade Kit as described in section 3.

2 Prerequisites



SPARE PARTS:

DX-D 100 Upgrade Kit, Spare Part Number CM+6056420, consisting of the following components:

- SparkLAN WUBR-507N WiFi network card
- DX-D 100 on/off board
- Installation instruction (WiFi on/off antenna)
- WIN7 upgrade license sticker on PC (WIN XP to WIN7)



NOTE:

The upgrade license is limited to a quantity of 40 pieces and can not be ordered anymore because only a limited number of XP systems (< 40) were delivered to the field.



SOFTWARE:

Download the T605220.0 Sparklan Driver RT2870_V3.2.1.0 from the HealthCare Library:

Direct Radiography → DR Equipment → DX-D 100 → Software



TOOLS:

Standard service engineers' tool kit

Flash storage device for storing the T605220.0 Sparklan Driver RT2870_V3.2.1.0



REFERENCED DOCUMENTS:

DX-D 100 - Clean Install Procedure, DD+DIS040.11E

Direct Radiography → DR Equipment → DX-D 100

3 Instructions

**REQUIRED TIME:**

Approximately 45 minutes for installing the Upgrade Kit, the driver for the SparkLAN WUBR-507N WiFi network card and the softconsole.

**CAUTION:**

Keep the unit turned off and isolated from the power supply.

3.1 Replacement procedure for the WiFi board kit:

- (1) Remove the frame that encloses the touch screen monitor.
- (2) Remove the rear cover.
- (3) Open the PC unit (pivot PC module on the right).

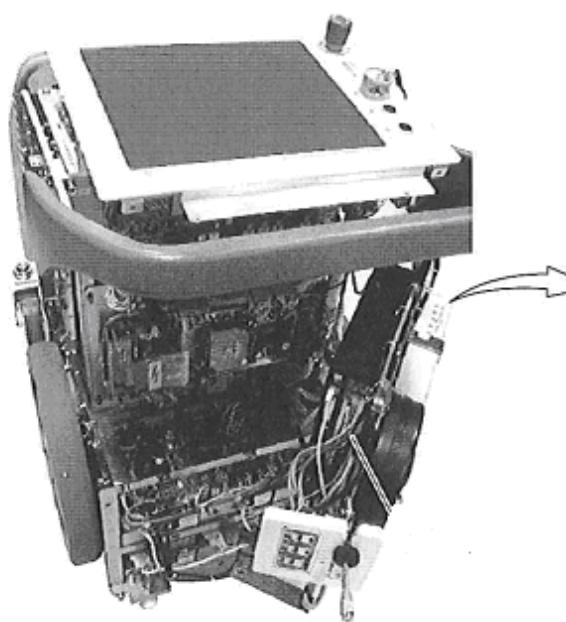


Figure 1

- (4) Once the back cover has been dismantled from the mobile unit, open the PC module and unplug all the PC connectors (see Figure 3 below).

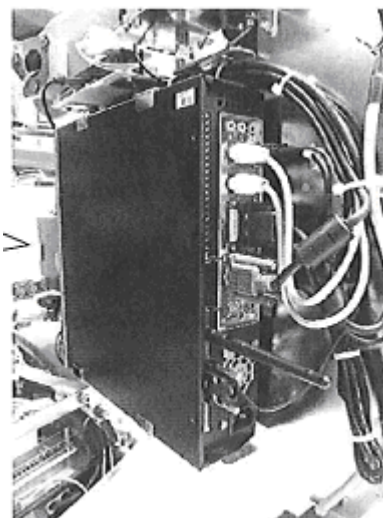


Figure 2

- (5) Loosen the screws of the PC fixing system to dismount the PC.
(6) Place the PC on a table and dismount the cover by removing the fixing screws (see figure 3).

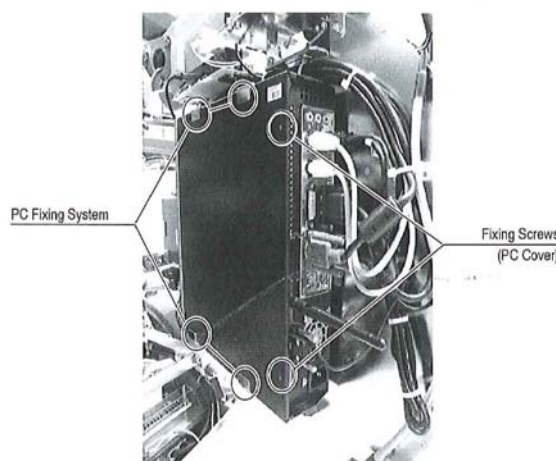


Figure 3

- (7) At the back of the PC, remove the fixing screw and dismount the fixing plate in order to remove the auto ON/OFF PCB in the next step.

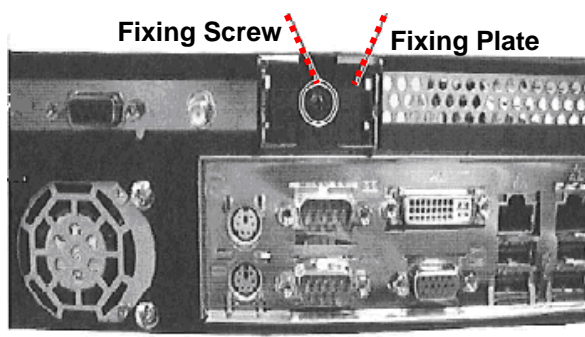


Figure 4

- (8) Inside the PC, remove the WiFi fixing screws, disconnect the power supply cables of the WiFi PCB and the auto ON/OFF PCB, then remove both boards.

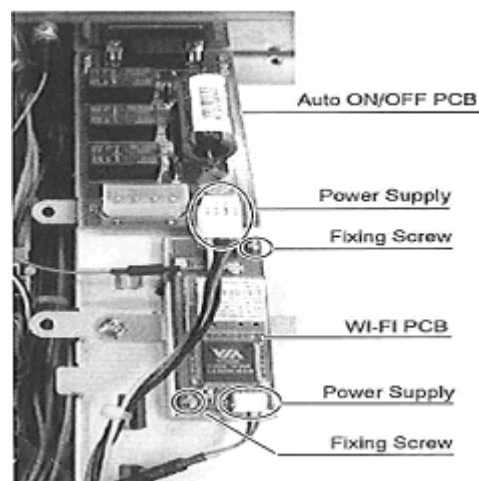


Figure 5

- (9) Install the PCB inside the PC.



Figure 6

- (10) Install the new auto ON/OFF mobile PCB inside the PC.
(11) Connect the power supply cables of each board.

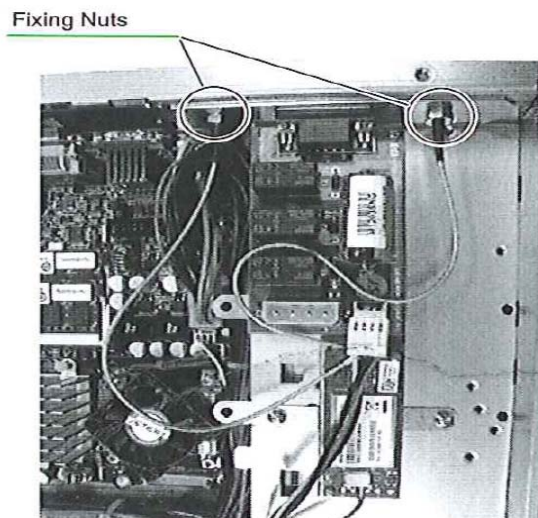


Figure 7

- (12) At the back of the PC housing, fasten the fixing plate.

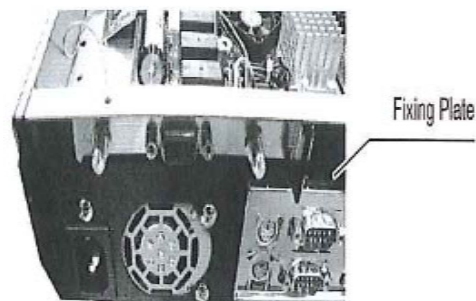


Figure 8

- (13) Mount the PC cover, install the PC in the PC module with the fixing system and connect all the pc connectors.
- (14) Close the PC housing.
- (15) Restart the PC.

After restart, install the driver for the SparkLAN WUBR-507N WiFi network card:

- (1) Insert the Flash drive medium at the PC.
- (2) Extract **Sparklan Driver RT2870.zip** from flash drive.
- (3) Double-click the exe inside: **IS_AP_STA_RT2870_D-3.2.1.0_VA-3.2.1.0_W7-3.2.1.0_RU-4.1.0.0_AU-4.1.0.0_051011_1.5.12.0WP_Free.exe**
- (4) Select **<install driver only>** and continue installation.

3.2 Installation of WINDOWS7

- (1) Install Windows 7 according to the DX-D 100 Clean Install Procedure (see 'Referenced Documents' above).
- (2) Put the Microsoft WIN7 license sticker on the PC when the upgrade has been performed. (License key is asked during install O.S).



IMPORTANT:

Check if the Sparklan Driver is still installed after finishing the Clean Install Procedure!

4 Verification

Check if performance and connection reliability of the DX-D100 are stable after installation of the Upgrade Kit.

Check if the IEEE802.11n protocol is supported after installation of the Upgrade Kit.

Check if the connection to the hospital network is working.

5 Keywords

WiFi, WLAN, Upgrade Kit, IEEE802.11n protocol

6 Version history

Version	Change	Date
1.0	Initial Version	08-2012



Details as of PDF Creation Date

Document Metadata

Title:	DX-D 100 Service Bulletin No. XX - New WLAN Adapter Available for XP-based DX-D 100
Livelink ID:	37440705
Version#:	10
Version Date:	2012/08/03 02:52 PM CET
Status:	Approved on 2012/08/06 03:51 PM CET
Owner:	Beate Richter (axnwp)
Created By:	Beate Richter (axnwp)
Created Date:	2012/06/25 06:40 PM CET
PDF Creation Date:	2012/08/06 03:51 PM CET

This document was approved by:

Signatures:

1. Josef Wagner (agwj) on 2012/08/06 03:30 PM CET
2. Paul Merckx (amdag) on 2012/08/06 01:40 PM CET
3. Bart Biesemans (amajm) on 2012/08/06 10:33 AM CET

Detailed Approver History:

- **Approval Workflow started on 2012/08/03 03:06 PM CET**
 - Approval task originally assigned to and completed by Paul Merckx (amdag) on 2012/08/06 01:40 PM CET
 - Approval task originally assigned to and completed by Josef Wagner (agwj) on 2012/08/06 03:30 PM CET
 - Approval task originally assigned to and completed by Bart Biesemans (amajm) on 2012/08/06 10:33 AM CET

Version & Status History

Version#	Date Created	Status
10	2012/08/03 02:52 PM CET	Approved - 2012/08/06
9	2012/08/03 02:51 PM CET	
8	2012/07/24 11:22 AM CET	Reviewed - 2012/07/31

7	2012/07/23 01:49 PM CET	
6	2012/07/23 01:28 PM CET	
5	2012/07/12 04:15 PM CET	
4	2012/06/27 06:01 PM CET	
3	2012/06/27 05:56 PM CET	
2	2012/06/25 06:49 PM CET	
1	2012/06/25 06:40 PM CET	