

Service Information Bulletin

This Bulletin is for information only.

Introduction of DMC firmware version 11R3b7

Task

Timing		Category	
Next service	<input type="radio"/>	Apply at all sites	
	<input type="radio"/>	Apply at affected sites as listed below	
	<input checked="" type="radio"/>	Optional to improve functionality of product	

Task Tracking

After completion of your task the following entry in your Service Report is required:

DD+DIS169.15E

*

* Insert the document number into the field "Comment" (SMS form).

Purpose of this document:

- This document contains information about updating to DMC Firmware V11R3b7.
- The new DMC firmware version improves the DX-D 100 driving behavior.

1 Introduction/purpose

The new DMC firmware version 11R3b7 improves the driving behavior of the DX-D 100 / DX-D 100 Wireless mobile unit at maximum speed.



IMPORTANT:

- In case a DMC board with revision “H” is already installed: Update the firmware to V11R3b7. Refer to section 3.1.
- In case a DMC board with a revision lower than “H” and defective firmware is installed or
- in case a DMC board with a revision lower than “H” and defective firmware is already installed: Order one of two available DMC board kits. Refer to section 0

2 Prerequisites



SOFTWARE:

Download the "DMC firmware V11R3b7" software from the Agfa HealthCare Library:
Direct Radiography > DX-D 100 > DX-D 100 Generic



SPARE PART:

For DX-D 100 / 100 Wireless units **with installed DMC board revision “H”**, order the following spare part:

- **SC+A3539-01** upgrade digital motion control V11R3b7, containing the DMC firmware 11R3b7.



NOTE:

SC+A3539-01 is an adaptation of the current spare part.

For DX-D 100 / 100 Wireless units with installed DMC board **with a revision lower than “H” and defective firmware**, order the following kit:

- **SC+A521079-03** Digital Motion Control V11R3b5, consisting of
 - DMC board rev H with firmware V11R3b5 already loaded
 - buzzer
 - instructions
 - tie wraps



IMPORTANT:

It is recommended to use this kit as it is and install the new DMC firmware version V11R3b7 on top. See sections 3.1 and 0 of this document.

For DX-D 100 / 100 Wireless units with installed DMC board **with revision “H” and defective firmware**, order the following kit:

- **SC+A521065-05** upgrade digital motion control V11R3b5 kit, consisting of
 - cable
 - CD with firmware V11R3b5, software and instruction

**IMPORTANT:**

It is recommended to use this kit as it is and install the new DMC firmware version V11R3b7 on top. See sections 3.1 and 0 of this document.



**NOTE:**

It is possible that spare part SC+A521065-05, which contains a DMC board with revision H and intact firmware, is in stock at some locations.

In this case it is possible to install this part. Additionally, it is necessary to upgrade the firmware version to V11R3b7.

- (1) Replace faulty parts with a spare part SC+A3539-01.
- (2) Upgrade the firmware. See section 3.1 of this document.

**REFERENCED DOCUMENTS:**

DX-D 100 - X-Ray Mobile - Troubleshooting, Document ID 41600560 ( [Agfa Intranet](#) /  [Agfa Portal via Internet](#)).

3 Instructions and verification

3.1 Upgrading the firmware with Digital Motion Control upgrade kit SC+A3539-01

**IMPORTANT:**

Upgrade of the firmware is necessary for DX-D 100 / 100 Wireless units **with installed DMC board revision “H”**.

**REQUIRED TIME:**

60 minutes

- (1) Upgrade the firmware version to V11R3b7 according to the description in: DX-D 100 - X-Ray Mobile - Troubleshooting - Chapter 6.8.6 Updating the Microcontroller, Document ID 41600560 ([Agfa Intranet](#) / [Agfa Portal via Internet](#)).

3.2 Replacing a DMC board with Digital Motion Control kit SC+A521079-03

**IMPORTANT:**

Replacement of the DMC Board is necessary for DX-D 100 / 100 Wireless units **with installed DMC board with a revision lower than revision “H” and for DX-D 100 / 100 Wireless units with revision “H” and defective firmware**.

**REQUIRED TIME:**

60 minutes

- (1) Replace the DMC board according to the description in DX-D 100 - X-Ray Mobile - Troubleshooting - Chapter 6.8.7, Replacement of the Digital Motion Control Board, Document ID 41600560, [Agfa Intranet](#) / [Agfa Portal via Internet](#).

**NOTE:**

Additional information for replacing the DMC board:

- Fix the DMC board with four plastic screws.
- Isolate the board from the metallic mounting plate by four plastic spacers.



Figure 1

- After replacement it is necessary to protect the DMC board SAT-A3539-01 against possible damages, e.g. during voltage measurement at Test Point TP8. For that, it is necessary to apply a coating protection over the C16 capacitor before measuring the voltage. Use Loctite 7400 for coating.
- In case Loctite 7400 is not available, use a piece of electrical tape over the C16 capacitor. In case the tape is used it should be removed after adjusting the TP8 voltage.

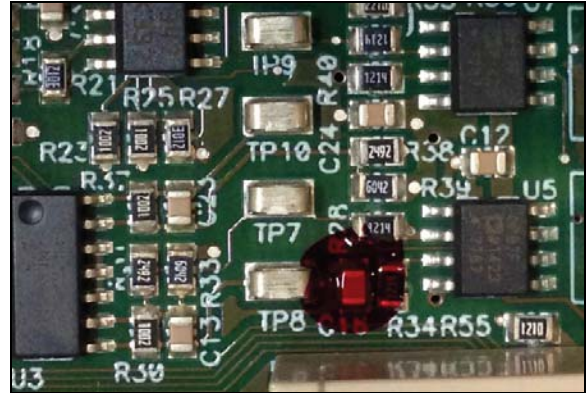




Figure 2

- (2) Upgrade the firmware version to V11R3b7 according to the description in: DX-D 100 - X-Ray Mobile - Troubleshooting - Chapter 6.8.6 Updating the Microcontroller, Document ID 41600560 ( [Agfa Intranet](#) /  [Agfa Portal via Internet](#)).

4 Keywords

DMC board, Firmware V11R3b7, revision "H"

5 Version history

Version	Change	Date
1.0	Initial version	11-2015