

Service Bulletin

Release of NX 2.0.8600 / 3.0.8600 and XRDI 10.1 for DX-D 100 System Version 1.0

Task			
Timing		Category	Scope
Next service as agreed with customer	<input type="radio"/>	Apply at all sites	PowerHelp complaint / HQ issue: HQ_1103310001
	<input checked="" type="radio"/>	Apply at affected sites as listed below	
	<input type="radio"/>	Optional to improve functionality of product	

Task Tracking

After completion of your task the following entry in your Service Report is required:

DD+DIS209.11E

*

* Insert the document number into the field "Comment" (SMS form).

Purpose of this document:

This document announces the introduction of the NX software version NX 2.0.8600 / 3.0.8600 and XRDI software version 10.1.

Affected serial number(s) / batch:

Installed DX-D 100 system 1.0 in combination with NX workstation version 2.0.8400 / 3.0.8400, XRDI version 10.0 and Varian software version L05.

The signatures on the approval page indicate the solutions described in this Service Bulletin have been reviewed and are NOT reportable because no actions are taken to reduce a "Risk to Health" according to our risk assessment process.

Use, dissemination, distribution or reproduction of this bulletin by unauthorized personnel is not permitted and may be unlawful. The controlled version of this document resides on Agfa HealthCare Library. Any printed copy of this document is uncontrolled.

Copyright © 2011 Agfa HealthCare N.V. | Document Node ID: 34250979 | service_bulletin_template_v06

1 Introduction

Symptom With NX software version 2.0.8300 / 3.0.8300 the following software issue has occurred:

Intermittent failure to display exposure parameters:

- If a patient is selected from RIS or manually entered the required routine is used for exam thumbnail to appear, i.e Chest PA.
- If the orientation is changed (Landscape/Portrait) from the original exam setting in the 'Edit' function within the 'Examination' window the exposure parameters do not appear within the GUI.

Solution Installation of the NX software version 2.0.8600 / 3.0.8600 in combination with XRFI software version 10.1. For more information refer to section 3 - Installation.

2 Prerequisites



SOFTWARE:

- Download the NX software version 2.0.8600 / 3.0.8600 from the FTP Server:



NOTE:

To access the FTP server, sign on using the username and password listed in the FTP Access document on the Agfa HealthCare Library, ([Intranet Link](#) / [Extranet Link](#)).

- For upgrades order the NX 2.0.8600 / 3.0.8600 upgrade kit (Starterkit DVD 1 and 2):
 - Via spare part code*: CM+6023930 (= Clean Installation Kit) **or**
 - Via Sales ABC code: 5QSNR

* The last digit in the spare part number indicates the spare part revision at release of this document. When ordering, the current revision of the spare part is delivered.

- Download the software **XRFI 10.1.005.zip** (File size approximately 22,8 MB) from the Agfa HealthCare Library to a portable storage medium:
<Direct Radiography → DR Equipment → DX-D 100 → Software>



TOOLS:

- An offline configuration tool for NX 2.0.8600 / 3.0.8600 is available on DVD 1 of the NX 2.0.8600 / 3.0.8600 software and on the FTP Server.
- Portable storage medium (checked to be virus-free, e.g. USB Memory Stick)

**REFERENCED DOCUMENTS:**

- DD+DIS056.10E / No. 01 - NX 2.0.8600 / 3.0.8600 - Release of NX 2.0.8600 / 3.0.8600 ([Intranet Link](#) / [Extranet Link](#))
- DD+DIS055.11E / Service Documentation, Chapter 7 Upgrade Installation - Installation Instruction for Upgrade from earlier NX Versions ([Intranet Link](#) / [Extranet Link](#))
- DR XREDI Service Manual, DD+DIS072.09E ([Intranet Link](#) / [Extranet Link](#))

3 Instructions

3.1 Installation of NX Workstation 2.0.8600 / 3.0.8600

**REQUIRED TIME:**

Approximately 5 hours

- (1) Install the NX software version 2.0.8600 / 3.0.8600.

For more information refer to the NX Service Documentation, Chapter 7 Upgrade Installation, as listed in section 2 - Referenced Documents.

After installation of the NX software version 2.0.8600 / 3.0.8600 install the XREDI software version 10.1; see section 3.2.

3.2 Installation of XREDI 10.1

**REQUIRED TIME:**

Approximately 20 minutes

**IMPORTANT:**

Do **not** remove the previous XREDI software version manually before starting with the upgrade installation.

- (1) Update the required license file (*.alf) otherwise the XREDI software cannot be installed.

**NOTE:**

For further information concerning obtaining or creating the ALF (Application License File), refer to the Licensing Service Manual DD+DIS012.06E ([Intranet Link](#) / [Extranet Link](#)).

**IMPORTANT:**

During installation of the licenses for XRDl software version 10.1 a warning will appear that that the current installed version does not match with the license file.

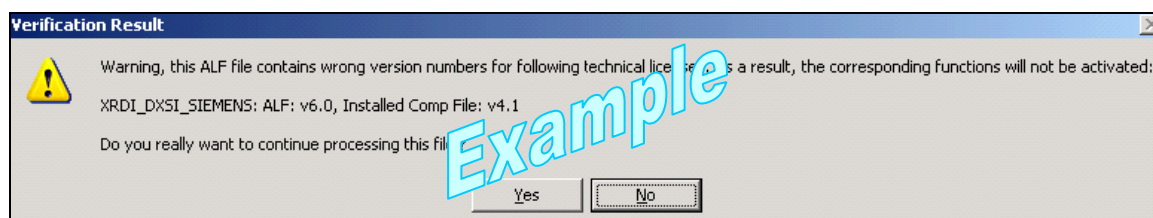


Figure 1

(2) Click **<Yes>** to continue the installation.

**NOTE:**

Following possibilities exist to install the XRDl software:

- For a manual installation of the XRDl software; perform step (3) to (6).
- For an alternative installation via NX; perform step (7) to (8).

(3) Stop the NX application before installing the XRDl Software version 10.1.

Go to: **Start → Agfa → NX → Stop NX**

(4) Install the XRDl software version 10.1 as described in the XRDl Service Manual; see referenced documents in section 2.

(5) After installation restart the NX application via:

Start → Agfa → NX → Restart NX Completely

(6) Go to **Start → Settings → Control Panel → Add or Remove Programs** and verify that the installed XRDl software version is: **10.1.005**

After a successful installation refer to section 4 - Verification.

(7) Put the XRDl **.exe** file in the folder **C:\ → Agfa → Healthcare → NX → Installers → XRDl** before running the NX upgrade.

**NOTE:**

During the NX workstation upgrade, the XRDl 10.1 software will be picked up by the upgrade process.

- (8) Make sure that the XRD1 software is selected as shown in Figure 2.

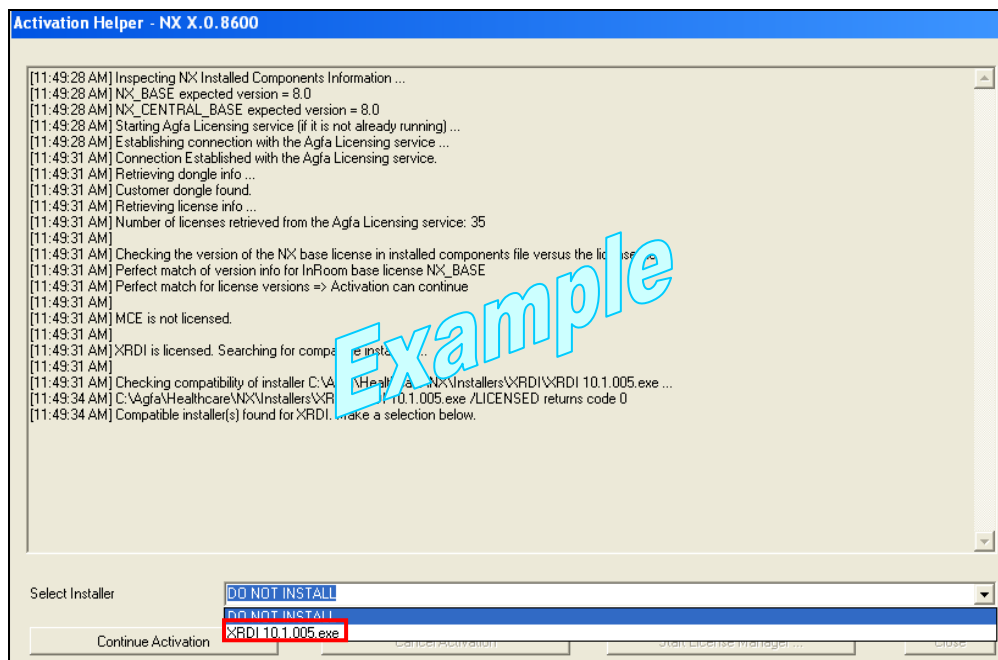


Figure 2

After a successful installation refer to section 4 - Verification.

4 Verification

The verification of the correct functionality of the installed DX-D 100 system in combination with NX Workstation version 2.0.8600 / 3.0.8600 is carried out by performing the following Acceptance Tests procedures:

- On system level: DX-D 100 System Acceptance Test
- On component level:
 - X-ray Modality Acceptance Test procedure
 - NX workstation Acceptance Test procedure
 - DR Detector

5 Keywords

- DX-D 100 system 1.0
- NX 2.0.8600 / 3.0.8600
- XREDI 10.1
- Upgrade Procedure

6 Version History

Version	Change	Date
1.0	Initial Version	09-2011



Details as of PDF Creation Date

Document Metadata

Title:	DX-D 100 - Service Bulletin No. 09 - Release of NX 2.0.8600 / 3.0.8600 and XRDI 10.1
Livelink ID:	34250979
Version#:	6
Version Date:	2011/09/07 08:22 AM CET
Status:	Approved on 2011/09/12 08:37 AM CET
Owner:	Claudia Tobi (axgzj)
Created By:	Claudia Tobi (axgzj)
Created Date:	2011/08/22 03:56 PM CET
PDF Creation Date:	2011/09/12 08:38 AM CET

This document was approved by:

Signatures:

1. Josef Wagner (agwj) on 2011/09/07 10:45 AM CET
2. Paul Merckx (amdag) on 2011/09/09 01:16 PM CET
3. Bart Biesemans (amajm) on 2011/09/07 08:50 AM CET

Detailed Approver History:

- **Approval Workflow started on 2011/09/07 08:26 AM CET**
 - Approval task originally assigned to and completed by Paul Merckx (amdag) on 2011/09/09 01:16 PM CET
 - Approval task originally assigned to and completed by Bart Biesemans (amajm) on 2011/09/07 08:50 AM CET
 - Approval task originally assigned to and completed by Josef Wagner (agwj) on 2011/09/07 10:45 AM CET

Version & Status History

Version#	Date Created	Status
6	2011/09/07 08:22 AM CET	Approved - 2011/09/12
5	2011/09/07 08:22 AM CET	
4	2011/09/02 10:34 AM CET	

3	2011/08/26 02:29 PM CET	Reviewed - 2011/09/02
2	2011/08/22 03:57 PM CET	Reviewed - 2011/08/26
1	2011/08/22 03:56 PM CET	